
arzum

SATEEN TOUCH

AR5176-02 IPL TUV ALMA SISTEMI
AR5176-02 IPL HAIR REMOVAL SYSTEM
AR5176-02 IPL HAARENTFERNUNGSSYSTEM
دليل النظام IPL AR5176-02 إزالة الشعر بالليزر

KULLANMA KILAVUZU
INSTRUCTION MANUAL
GEBRAUCHSANWEISUNG
دليل المستخدم

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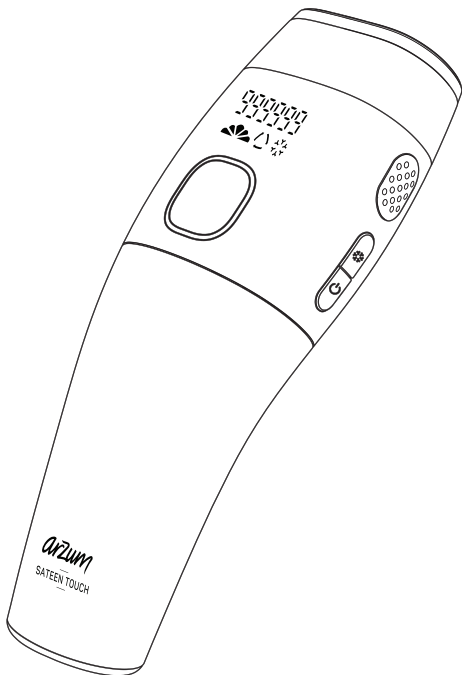
0850 222 1800

www.arzum.com

AR5176-02 ARZUM SATEEN TOUCH IPL HAIR REMOVAL SYSTEM

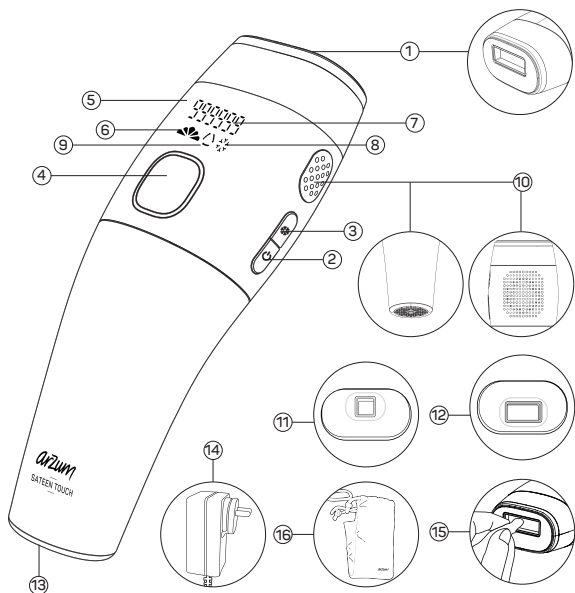- 1- UV-filtered light window
- 2- On/Off and energy level adjustment button
- 3- Cooling function button
- 4- Light pulse button and LED indicator
- 5- LED display
- 6- Energy level indicator
- 7- Remaining light pulse indicator

- 8- Cooling function indicator
- 9- Auto mode indicator
- 10- Air channels
- 11- Magnetic face cap
- 12- Magnetic bikini cap
- 13- Adapter input
- 14- Adapter
- 15- Cleaning cloth
- 16- Travel bag

Thank you for purchasing **Arzum IPL Hair Removal System**. This is one of the products designed to make life easier for you.

Before using the device, carefully read the instructions for use. Keep the user manual, warranty certificate, receipt, or invoice, and if possible, the box as well. If you give this device to someone else, please give them the user manual along with it.

SAFETY INSTRUCTIONS

- Use the device only for removing unwanted hair. After prolonged use, a reduction in the number of growing hairs may be observed. Use for any other purpose is inappropriate and may be dangerous.

- Do not use on animals.

- Do not use the device on any surface other than the skin.

- Do not use the device for commercial purposes.

- When not in use, during cleaning, or in case of malfunction, turn off the device and unplug it.

- If you will not be using the device again, it is recommended that you render it unusable by removing all dangerous parts, especially for the safety of children who may be tempted to play with the product.

- Provided that training is given on safe use and potential hazards, the device can be used by children aged 8 and above, as well as individuals with physical, sensory, or mental disabilities, or those lacking experience and knowledge. However, children should not play with the device; they should not be allowed to clean or maintain the device on their own.

- If your device is damaged for any reason, do not use it, do not open it, and do not attempt to repair it yourself. In case of any malfunction, take

it to an **Arzum Authorized Service Center**.

- For safe and efficient operation of the device, use only original parts and adapters.
- Do not use products with damaged cables. Contact an Arzum Authorized Service Center.
- Do not allow children to play with the product packaging. Plastic bags should not be left where children can reach them due to the risk of suffocation.
- Before connecting the device, ensure that the main voltage information matches the voltage information specified in the product's technical specifications. The technical specifications are listed on the device and on the adapter.
- Protect the device from sunlight, rain, and moisture.
- Do not use the device near water or on wet surfaces (bathtub, shower, sink, etc.).
- Do not use the device if it gets wet. Using a wet device may cause electric shock.
- Do not use the device if the light window (1) is cracked or missing.
- Do not touch the light window during or immediately after use; it may be hot.
- Do not use petroleum-based or flammable cleaning agents (alcohol, acetone, etc.); they pose a fire hazard.
- Do not cover the cooling vents. The device may overheat, posing a risk of electric shock or fire.

General health status and chronic illnesses

Do not use the epilator device if you have any of the following health conditions, or consult your doctor beforehand:

- If you are under 18 or over 65 years of age.
- If you have diabetes or any other systemic or metabolic disease.
- If you have congestive heart disease.
- If you have active implants such as a pacemaker, neurostimulator, or insulin pump.
- If you have had a disease that suppresses the immune system (e.g., HIV infection or AIDS).
- If you are taking medications that suppress the immune system.

- If you have any clotting disorder.
- If you have photosensitive epilepsy.

Sensitivity and Skin Sensitivity Conditions

- If you have a skin tone that is not suitable according to the "Suitable skin tones" table in the user manual.
- If you have a condition related to light sensitivity, such as polymorphic light eruption (PMLE), solar urticaria, or porphyria.
- If your skin is sensitive to light and easily develops rashes or allergic reactions.
- If you are using medications or substances that increase light sensitivity and their package inserts indicate a risk of photoallergic or phototoxic reactions.
- If you are using painkillers that reduce the skin's sensitivity to heat.

Medical History

- If you are pregnant or breastfeeding.
- If you have used Accutane, Roaccutane, or any other form of isotretinoin within the last 6 months.
- If you are undergoing skin treatment with Alpha Hydroxy Acid (AHA), Beta Hydroxy Acid (BHA), topical tretinoin, or azelaic acid, or if such treatment was administered within the last 7 days.
- If you are taking blood-thinning medications, including high doses of aspirin, and require at least a 1-week washout period for these medications.
- If you have undergone radiation therapy or chemotherapy within the last 3 months.

Skin Conditions and History

- If you have a history of collagen disorders, such as keloid formation or slow wound healing.
- If you have active skin cancer, a history of skin cancer, or a history of local cancer in the areas to be treated.
- If you have a skin pigmentation disorder, such as vitiligo.

Conditions Related to the Treatment Area

- If you have a history of vascular disease in the treatment area, such as varicose veins or vascular dilation (vascular ectasia).
- If there is an infection, eczema, burn, folliculitis, open wound, abrasion,

cold sore (herpes simplex), lesion, or bruise (hematoma) in the areas where the application will be performed.

- If the skin is irritated (redness, cuts), sunburned, tanned, or if a tanning product has been applied.
- On moles, freckles, large veins, dark pigmented areas, scars, or skin imperfections without consulting your doctor.
- On areas with warts, tattoos, or permanent makeup.
- Around the eyes and eyebrows. It can only be used on the body and on the female face (under the cheekbones).
- On the lips, scalp, nipples, between the legs, around the nipples, labia minora, vagina, anus, nostrils, and inside the ears.
- Areas where long-lasting deodorant is used. This may cause skin reactions.
- On or near artificial objects such as silicone implants, subcutaneous injection channels (e.g., insulin dispensers), or piercings.

For men: Shoulders and below (chest, back, arms, abdomen, and legs)

WARNING: Do not use on the face, neck, or genital area. Use on facial hair or beard may cause permanent or uneven results; this may lead to changes in physical characteristics.

NOTE: This list does not cover all possibilities. If you are unsure whether you can use the device, be sure to consult your doctor.

IS IPL HAIR REMOVAL TECHNOLOGY SUITABLE FOR YOU?

For IPL (Intense Pulsed Light) technology to be effective, there must be a clear contrast between skin tone and hair color. The best results are achieved by users with fair skin and dark hair. IPL light is absorbed by melanin, the natural pigment found in skin and hair. This absorption heats the hair follicle, rendering it ineffective.

However, effectiveness may decrease or risks may arise in the following cases:

- **Red, white, gray, and light yellow hair** may not absorb enough light due to their low melanin content, and the desired results may not be achieved.
- **Very dark skin tones** may absorb too much light due to high melanin content, which can cause burns on the skin.

NOTE: For first-time users, energy level 1 is recommended to assess skin sensitivity. After the skin test is completed, wait 24 hours and check for any reactions. Higher energy levels provide more effective results. For subsequent applications, it is recommended to use the device at an increased energy level.

SKIN TONE & HAIR COLOR CHART TREATMENT SCHEDULE

Please check your skin tone and hair color against the table before using the device to ensure compatibility. Using the device with incompatible combinations may cause skin irritation, burns, or ineffective results.

TEN RENGİ						
TÜY RENGİ						

APPLICATION SCHEDULE



Recommended frequency of use for faster and more effective results:

- **Weeks 1–2:** Apply 3 times per week.
- **Weeks 3–4:** Continue applying 2 times per week.
- **Weeks 5–8:** Apply 1–2 times per month.

NOTE: After week 8, if additional hair growth is observed, maintenance is recommended every 1–3 months.

The time between treatments varies depending on hair growth rate and body area; effects may vary from person to person.

Visible Results After 8 Weeks:

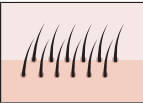Week 1		You will still see some hairs growing. These are likely to be hairs that were missed during treatment.
Week 2		You should see a reduction in hair growth. However, many hairs may still not have been treated in their growth phase.
Week 4		You should see a significant reduction in hairs within the treated area.
Week 8		Any remaining hairs should be finer and lighter in color.

MAGNETIC HEADS

- To achieve the best results, it is important to use the correct head designed for the specific body area you want to treat.
- In addition, thanks to their magnetic structure, the heads can be easily attached and removed, providing practicality and comfort during use.

This product includes two magnetic heads:

Bikini Area Head

- Provides effective and safe application on the bikini area, which has a more sensitive and dense hair structure. It offers comfortable use by establishing more controlled contact with the skin.


Face Area Head

- Specially designed for sensitive areas such as the face. Its small application surface allows for precise and safe light emission on areas like the upper lip, chin, and cheeks.
- Since these heads are used on different body areas, it is important to use them separately for hygiene reasons.

USE OF THE DEVICE

Preparation Before Application

- Prepare the skin by shaving the hair in the area to be treated, then clean and dry the area.
- This step allows the IPL light to reach the hair follicles directly and increases the effectiveness of the application. Methods other than shaving may remove the hair follicle, which can negatively affect the performance of the application. If hair remains on the surface, a burning smell will occur during the application. Therefore, it is not recommended to remove hair with gel, cream, or an epilator.

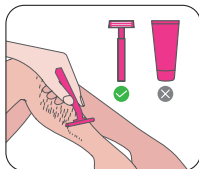Operating the device

- First plug the device's adapter into the power outlet, then plug the adapter socket into the device. The device only operates when connected to the adapter. When connected, "OFF" appears on the display.

- Turn on the device by pressing and holding the On/Off button for 2 seconds. The fan will start running. The device may heat up slightly during use.
- The device has 5 different energy levels. Press the On/Off button (2) to adjust the energy level. When the energy level is selected, the icon indicating the level will appear on the screen. (🌞)

Cooling function: It provides a more comfortable experience during application by cooling the skin. To activate this feature, press the cooling function button (3). The product has two different modes: automatic and manual.

Manual Mode: Press the light pulse button (4) briefly to perform a single light pulse each time. Suitable for small and sensitive areas such as fingers and upper lip.


Automatic mode: To activate automatic mode, press and hold the light pulse button (4) for 2 seconds. The automatic mode icon appears on the screen. (⏏)

- To start the procedure, place the device vertically on the skin. Then, slide the light window over the skin; the device automatically emits light when it detects contact. This way, you don't need to press the

light emission button repeatedly. This mode is ideal for large application areas such as arms and legs. To turn off the function, press and hold the light pulse button for 2 seconds; the icon on the screen will disappear.

- You can continue the procedure by selecting one of the above modes.

- Press the light output window (1) so that it is in full contact with the area to be treated. When the window is in full contact with the skin, the LED around the light pulse button will start flashing; this indicates that the device is ready to emit light pulses. Press the button to emit a light pulse. A slight warmth will be felt on the skin. Move the light window to the next area and repeat the process.

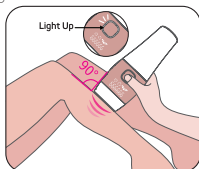- Ensure that light pulses are applied to the entire area, but do not apply light pulses to the same area twice (avoid overlapping light pulses).


- When the process is complete, press and hold the power button. The fan will continue to run for approximately 10 seconds to cool the device; during this time, a countdown timer will be displayed on the screen, and then the device will automatically shut down.

WARNING: The device will automatically shut down if not used for 3 minutes.

When the internal temperature of the device reaches 75°C, an error message will appear on the screen.


Please cool the device in a well-ventilated area, after which it will resume normal operation.

- Unplug the device from the outlet and power source.

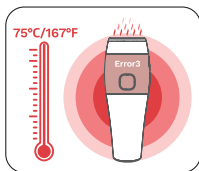CLEANING

- Before cleaning the device, make sure it is turned off and unplugged.
- The light window may become hot during use; ensure it has cooled completely before cleaning.
- Do not use dish sponges, abrasive cleaning materials, gasoline, acetone, or other liquids that may damage the surface.

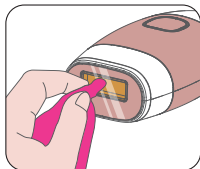- After use, wipe the light window with the special cleaning cloth provided with the product
- If necessary, wipe the exterior of the device with a clean, soft cloth.
- To get the best performance from the device, clean it after each use.

SIMPLE TROUBLE-SHOOTING

When fault occurs during the usage of the instrument, users can follow the following guidance to eliminate some simple faults. If the problem can not be solved after all steps are finished, please contact the Customer Service Department so as to solve faults as soon as possible.

Problem	Solution
The device does not power on when the Power button is pressed.	Verify the device's electrical connection. Ensure that the adapter is securely connected to the device socket. Confirm that the adapter is plugged into a live electrical outlet.
The device does not emit a light pulse when the Light Pulse button is pressed.	· Check whether the device has overheated. The device will not emit pulses if it is too hot. · Ensure that the treatment window is in full contact with the skin and that the Light Pulse READY indicators on both sides of the device are illuminated.
The energy level does not change when the Energy Level button is pressed.	Pressing the Energy Level button repeatedly or rapidly may trigger the device's 3-minute safety lock. Wait for 3 minutes, then press the button again to select the desired energy level.
Significant discomfort is experienced on the skin during device use.	• Hair trimming: Shave the treatment area prior to use to shorten hair. • Recent tanning: Exposure to sunlight or artificial tanning products shortly before use may increase skin sensitivity. Avoid sun exposure and tanning devices for at least 2 weeks prior to treatment. • Energy level adjustment: Reduce the energy level to a comfortable setting.
The device is powered on, lights are illuminated, but no pulse is emitted.	The device's safety system may have activated due to overheating. Allow the device to cool; normal pulsing will resume. This is a safety feature and does not indicate a malfunction.

Model No	AR5176
Light Output Area (Spot Size)	3 cm ²
Wavelength	600–1200 nm
Total Energy	15J
Energy Density	5 Levels
Repetition Rate	1.4–2.8 seconds (depending on the selected energy level)
Operating Mode	Auto Mode / Manual Mode
Product Size	19.51 × 7.23 × 4.763 cm
Additional Header (Spot Size)	1 cm ² / 2 cm ²

Product Life 7 years

Technical Specifications

Main Unit: 12.0 V $\equiv$ 4A

Adapter Model: DZ048BHL120400V

Adapter Input: 100-240V~50/60Hz 1.5A

Adapter Output: 12.0V $\equiv$ 4.0A 48.0W

HANDLING AND SHIPMENT

- During handling and shipment, your appliance must be kept in original package in order to prevent damages to its parts.
- Keep at normal position during shipment.
- Do not drop the product when transporting and protect from impacts.
- After delivery of the product, failures and damages that may occur while transporting are not covered by warranty.

The product is in compliance with the WEEE Directives.

WEEE

This product contains recyclable materials that compliance with the WEEE directives. Do not dispose this product as unsorted municipal waste. Please contact your local municipality for the nearest collection point.

Manufacturer:

Foshan Jindi Electric Appliance Co., Ltd

No. 13 Baiye Avenue, Xiqiao Science and Technology Industrial Park, Xiqiao Town, Nanhai District, Foshan City, Guangdong Province, China
0086 75786868180

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