# arzum

## SMARTFIT NEO


AR 5090 AKILLI VÜCUT ANALİZ BASKÜLÜ  
AR 5090 SMART BODY ANALYSIS SCALE  
AR 5090 SMART KÖRPERANALYZEWAAGE  
میزان تحلیل مکنونات الجسم الذكية AR 5090

**KULLANMA KILAVUZU**  
**INSTRUCTION MANUAL**  
**GEBRAUCHSANWEISUNG**  
دليل المستخدم

**3**  
**YIL**  
**YEARS**  
**JAHRE**  
GARANTI  
WARRANTY  
GARANTIE

# SMARTFIT NEO AR 5090 AKILLI VÜCUT ANALİZ BASKÜLÜ

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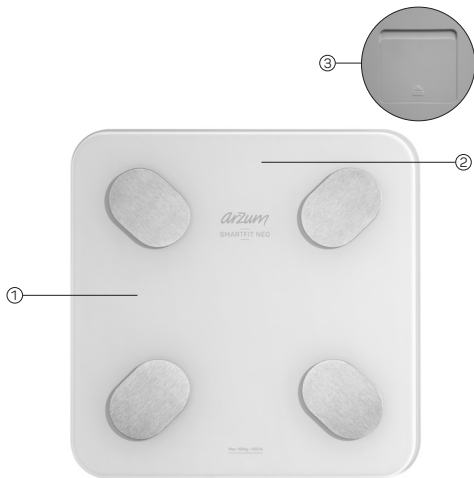
**www.arzum.com**

# SMARTFIT NEO AR5090 SMART BODY ANALYSIS SCALE

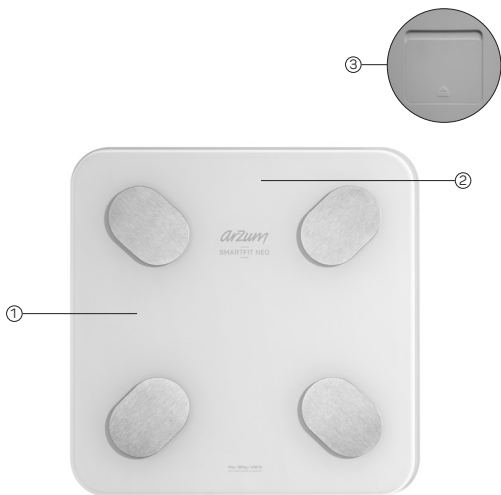- 1- Ultra transparent glass surface
- 2- LED display and comfortable readability
- 3- Battery section

## Features

- Extremely precise sensor system
- Up to 180kg/400 lb
- 50 gr/0.1 lb accuracy
- Weight unit kg, lb
- Changing kg, lb over Wireless connection
- Thanks to the wireless connection application; BMI analysis, weight tracking, body size tracking, water ratio tracking can be done.
- Automatic shutdown
- Automatic reset
- Low power indicator
- Overload indicator

We would like to thank you for choosing **Arzum scale**. We hope that this product will make your life easier like other Arzum products.

We recommend you to read this instruction manual carefully before using your scale and keep the manual as it may be necessary in the future.

**ARZUM SMARTFIT NEO SMART BODY ANALYSIS SCALE** is designed for domestic use. Do not use the appliance for commercial purposes and trading.

### Warning:

- Data of evaluating the body fat and hydration in percentages are only the reference values. It is recommended not to use it for any medical purpose.
- It is recommended for persons that use heart pacer or other medical equipments having electronic components do not use this device.
- Do not use your device in the places where there is some equipment which could produce high electromagnetic waves in order to prevent attraction. Use your scale either at an appropriate distance or by shutting down such equipments.


## Measurement Tips

- This product uses BIA (bioelectrical impedance analysis) technology. In order to obtain correct and consistent measurement results, please

use it at the same time during the day and on the same conditions as possible. For accurate results, do not measure in the following situations.

- During or after exhausting exercises
- Immediately after taking a bath or after steam bath
- During alcohol consumption and fever
- Stand upright on bare feet during measurement. Measurement will not be accurate if socks are worn. Dirty feet can also affect measurement accuracy.
- The weight unit may be set as pound (lb.) or kilogram (kg) with the mobile application. The height unit will also change accordingly (meters vs. foot/inch).

Place the scale on a hard floor to measure the weight correctly.

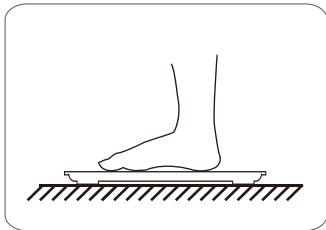- Accurately measure body fat and other body composition parameters: Stand upright on the scale and make sure that both feet touch the electrodes equally, as shown below.

## Recommendations of Usage and Maintenance

- It is recommended to use **Arzum Smartfit Neo smart body analysis scale** on a hard and flat surface.
- Always measure your weight with the scale placed on the same ground. Uneven ground can affect weight figures.
- Measure your weight before meals and always at the same time, without your clothes and without shoes. To get your body analysis

values step on the scale without socks.

- Clean the scale surface with a damp cloth, but do not allow water to get inside. Do not use chemical cleaning agents.
- Do not store your scale in highly humid places, do not put it in water or other liquids.
- Do not store your device in an upright position when you are not using it.
- Do not drop your scale and do not jump on it.
- In case of an issue that you cannot solve, please contact your nearest Arzum Authorized Service.
- Use your device as recommended and fit for the purpose.

## Battery Installation

Install three 1,5 volt AAA batteries in the battery section under the scale.


## OPERATING YOUR DEVICE

### Normal weighing:

- Ensure that the scale is placed on a firm and flat surface.
- Step on the scale so that your feet are evenly spaced and wait until the weight is measured without moving.
- The value measured will be displayed on the screen thanks to the sensitive sensors in a short time.

### Overload indicator:


When the scale is overloaded (180 kg and more) the display will show "EEr".
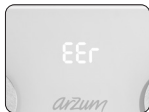
### Low power indicator:

Smartfit scale will give a "Lo" warning on the display when the battery level reaches a low level.

If you see the "No Display" warning, please change the battery of your scale.
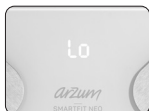
## WEIGHING WITH WIRELESS CONNECTION

### Download Application:

You can download the **FeelFit app** on your tablet or mobile phone from "App Store" or "Google play".


### Registration:

After downloading the application, you will be asked to register on the screen to be opened.

Enter your e-mail address to create your password.


Complete the other steps in order.
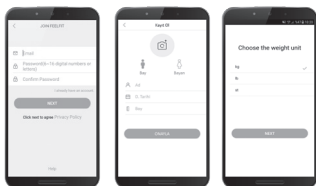
### Connect the device:

Turn on the wireless feature of your smart device.

· Open the FeelFit application go to, wireless Connections on your smartphone and find your scales.

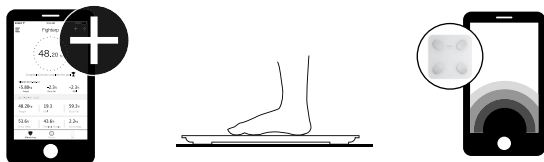STAND UPRIGHT ON BARE FEET DURING MEASUREMENT. DEVICE WILL NOT PROVIDE CORRECT RESULTS IF SOCKS ARE WORN.

## DISPLAYING BODY ANALYSIS

### Analysis report :

You can easily access the following details in the body analysis.

- **Weight:** One of the important signs of personal health.
- **BMI:** Shows the body mass index.
- **Body fat percentage:** Shows the percentage of fat in the body.
- **Body type:** General description and evaluation of body shape of the person.
- **Fat-free body weight:** One of muscle weights except for fat.
- **Subcutaneous fat percentage:** The percentage of fat tissue stored in the skin to body weight.
- **Visceral fat accumulation:** Fat accumulation in the body different from subcutaneous fat; It occurs mainly in the abdominal cavity, occurring around organs in humans.
- **Body Water percentage:** Percentage of water in human body including blood, lymph, extracellular fluid, etc. in terms of weight.
- **Skeleton muscle:** The ratio of muscle involved in the mechanical system of our limbs and other parts of the body.
- **Muscle mass:** The total muscle weight, including skeletal muscle, cardiac, and smooth muscle.
- **Muscle storage ability level:** The size of muscle storage capacity represents the number of human muscles and reflects the current level of the body's ability to retain muscle.
- **Bone mass:** Bone tissue consists of bone minerals and bone matrix per unit volume.
- **Protein:** Plays a vital role in the body. Builds and protects muscles, organs and other tissues.
- **Basal metabolic rate:** The least level of energy needed to survive in an inactive state.
- **Metabolic age:** Ideal physical body age is 2/3 of the actual age.

## Measurement:

1-Step on the scale and wait for the value to stabilize.  
Compare the weight value from your smart device.

2-When zero figures on the screen flash and the last value is stabilized when you last stepped on the scale, the body analysis will be completed. You can check this from your smart device.

## Adding New Users:

You can create a new record by going to "Add Member", "Add Measurement User" in "My Account" menu in FeelFit application.

If there are more than one users, you should select correct member manually. Otherwise, app will warn if there is big difference than previous measurement.

YOU CAN CHANGE APP LANGUAGE TO YOUR PREFERRED LANGUAGE THROUGH SETTINGS AT ANY TIME.

Product life 7 years  
(3X1,5V AAA)

## HANDLING AND SHIPMENT

- During handling and shipment, your appliance must be kept in original package in order to prevent damages to its parts.
  - Keep at normal position during shipment.
- Do not drop the product when transporting and protect from impacts.
- After delivery of the product, failures and damages that may occur while transporting are not covered by warranty.

**The product is in compliance with the EEE Directives.**

## WEEE

This product contains recyclable materials that compliance with the WEEE directives. Do not dispose this product as unsorted municipal waste. Please contact your local municipality for the nearest collection point.

**Manufacturer:**

**SHENZHEN YOLANDA TECHNOLOGY CO., LTD.**

2F JINFULAI BUILDING NO.49-1, DABAO ROAD, BAOAN 28 DISTRICT, SHENZHEN

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