# arzum

**SATEEN  
TOUCH PRO**

AR 5050 IPL TUV ALMA SİSTEMİ  
AR 5050 IPL HAIR REMOVAL SYSTEM

**KULLANMA KILAVUZU  
INSTRUCTION MANUAL**

**3**  
**YIL**  
**YEARS**  
GARANTİ  
WARRANTY

**Turkey**

# ARZUM SATEEN TOUCH PRO AR 5050 IPL TÜY ALMA SİSTEMİ

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Made in China-Menşei Çin'dir. Tel: 0850 222 1 800

**www.arzum.com**

# ARZUM SATEEN TOUCH PRO AR5050 IPL HAIR REMOVAL SYSTEM

- 1- Light Exit Window with Integrated UV Filter
- 2- Skin Contact Surface with Integrated Safety System
- 3 – Intensity Lights
- 4 – Power Light
- 5 – On-off/POWER Button

- 6- Adaptor Input
- 7 – READY Light
- 8 – Adaptor
- 9- Skin Tone Test Sensor
- 10-Flash Button
- 11- Treatment Glasses with UV Filter
- 12- Velvet Pouch

Thank you for purchasing **Arzum IPL Hair Removal System**. This is one of the products designed to make life easier for you.

## **SAFETY INSTRUCTIONS**

Read the operating instructions carefully before putting the appliance into operation and keep the instructions including the warranty, the receipt and if possible the box with internal packaging. If you give this device to other people please also pass on the operating instructions.

- This appliance must be used only for the purpose for which it was designed, namely as a treatment for hair removal for home use. Any other use is considered not compliant and therefore dangerous.
- The device is not intended for commercial use.
- Always disconnect the appliance from the mains supply before cleaning or performing maintenance and when the appliance is not in use or malfunctions.
- If you decide to stop using the appliance, it is advisable to make it inoperative by removing all potentially dangerous parts, especially to protect.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

- If the appliance is faulty or malfunctions, switch it off and do not tamper with it. Do not repair the device by yourself. Please contact an authorized technical service center.
- Do only use original parts.
- If the power cable is damaged, don't use the product and contact an authorized technical service center, so as to prevent all risks.
- Do only use original adaptor.
- Do not let children play with the packaging! Keep the plastic bag out of the reach of children – risk of suffocation!
- Before connecting the appliance, please make sure that the voltage data indicated on the identification technical data correspond to those of the mains voltage. The identification technical data are found on the appliance and on the battery charger.
- Do NOT expose the appliance to humidity or weather conditions (rain, sun, etc.)
- The device must only be used on human hair. Do not use the device on animals.

## **CONTRAINDICATIONS**

- Not for pregnant women, or in menstrual period, in suckling period or people with scar diathesis,
- Not for patient with heart disease, hypertension, malignancy, infectious skin disease, skin with burns and anti-coagulopathy,
- Not for person with photosensitization history or diabetes,
- Not for person who has depilated in other methods in the past 6 weeks,
- Not for person used skin whitening products or products with vitamin A in the past 1 week,
- Not for person who has taken essential oil treatment at the same time,
- Not for person who has taken Aspirin, immunosuppressant and used photoaesthesia drugs (like tetracycline, sulphonamide),
- Not for person who has been exposed to strong sunlight in the past 2 weeks,
- Not for sensitive areas, like mucosa etc, and private positions, like upper eye-lid, lip, nipple and vagina etc,

- This appliance is not intended for children under age of 12 years.

If you take prescribed medication or receive regular medical monitoring because of a certain condition, consult your doctor before you use Arzum SATEEN TOUCH PRO. You should also consult your doctor if you are not sure if the appliance is suitable for you.

## **INSTRUCTIONS FOR SAFE**

- Wear glasses during the operation. Do not look at the twinkling light or light to eyes directly.
- During the treatment, it is normal to have some skin redness, suggest to use ice pad after treatment.
- Before operation, use alcohol cotton pad or swab to disinfect the treatment window, but do not spray alcohol directly on the treatment window.
- To avoid destroying the device, please do not apply gel or other cosmetics on treated areas during the treatment.
- Try not eating light-sensitive food, like spinach, cole, mustard etc. during the treatment.
- Adjust the energy level from low to high. On condition of safety, the higher energy level, the better result. For dark skin, the skin reaction may be delayed. Thus shall not raise the energy level too quickly.
- Use skin moisture and sun block product after phototherapy treatment.
- Avoid exposing by ultraviolet light in using period, protect skin from the sun and avoid the computer radiating with skin moisture and sun block product in daily work.
- During the use, do not shield the air intake and louver to impact heat dissipation.
- In the operation, if the hair on bolt pastes on the plastic of the lighting window edge, will cause the plastic of the lighting window edge black or black spot, but cannot affect product function and normal use.
- In the latter period of specified life, it's normal that the in wall of lamp turns black, it's the aging process for lamp. The weaken rate of power is still in the specified scope, the device still can work normally.
- The product is not water-proof, so please do not use in bathroom, shower or any other places in which the instrument may be watered.

- In case that the instrument is soaked please power off before touching the instrument for your safety.
- Avoid keeping it in places where there is a strong magnetic field or electric wave occurred, or places in high temperature.
- Only use on skin, not on any other surface.

## LIFE TIME REMINDER SYSTEM

- There is a built-in life time reminder system in device and the total life time is 300.000 flashes.
- When shots emitted more than 50,000 shots (including 50,000 shots), the indicators of two sides will flash once.
- When turn on the device each time; when shots emitted more than 100,000 shots (including 100,000 shots), the indicators of two sides will flash twice when turn on the device each time.
- When shots emitted more than 150,000 shots (including 150,000 shots), the indicators of two sides will flash 3 times when turn on the device each time.
- When shots emitted more than 200,000 shots (including 200,000 shots), the indicators of two sides will flash 4 times when turn on the device each time.

## SKIN TONE TEST SENSOR

Function: Skin tone test sensor is used to test skin tone type and offers power level suggestion.

### Instruction:

- a) Voltage: 3 V d.c; Power; 0,3W
- b) Battery inside, no need charging,
- c) Use cotton swab or cleaning cloth to clean and disinfect the test window,
- d) Let the test window touch treated areas or the areas with the similar skin color flatly then press the turn on/off button,
- e) The skin type display will show what skin type by indicator,
- f) Skin tone test for more than 2 times to make sure test results accurately,

The skin tone sensor is only applicable for only **Arzum SATEEN TOUCH PRO**. Once you test your skin type, you can choose suitable energy level from Treatment Level and Times Comparisons Table for the treatment area.

**NOT:** Over temperature protection systems automatically becomes activated in case of overheating and device stops pulsing for a while.

This is not because of any defect on device, but just because of over temperature protection system of the device.

## PATCH TEST

- 24 hours before treatment check for any reaction by treating a small area first (Consult the Treatment Level and Times Comparisons Table),
- For initial operation, choose a small area on arm, treat the area in each power level up to the maximum power level for your skin type (repeat 1-2 shots on the same area),
- If there is any redness and swollen, burnt scares, change of hair follicles, irritation or any other discomfort, do not use the instrument before a thorough consultation with doctors,
- If no skin abnormality after 24 hours, feel strong hotness, then the last low level is the most proper energy for your skin. Please do regular treatments.

## PRODUCT INTRODUCTION

### Packing

Hair Removal, Skine Tone Test System, Glasses, Adapter, User Manual, Velvet Pouch

### Main Structure

Light Exit Window with Integrated UV Filter  
Flash Button  
Intensity Lights  
Power Light  
On-off/POWER Button  
Adaptor Input  
READY Light  
Adaptor  
Skin Tone Test Sensor

## PRODUCT FUNCTIONS INSTRUCTIONS

Application: Lip hair, armpit hair, body hair and legs, hair on places that affecting looks such as hair line on forehead and bikini area.

For the group with black or dark brown hair, not for white, grey, red or light blonde hair.


Arzum SATEEN TOUCH PRO can be used for 5 of 6 skin types below.

### Skin Color Comparasion Table

| I | II | III | IV | V | VI | Legend                                    |
|---|----|-----|----|---|----|---|
| ✓ | ✓  | ✓   | ✓  | ! | X  | ✓ Suitable<br>! Caution<br>X Not suitable |

## OPERATION STEPS

- 1) Use cotton pad or swap to disincect light exit window,
- 2) Connect to the power,
- 3) Press the power button for a while to turn on the device,

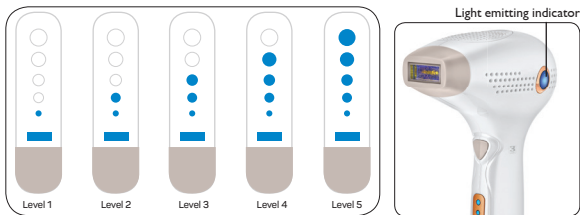4) Press power adjusting button to choose suitable power level,

**NOTE: Adjust from level 1 to level 5; level 1 is the weakest, level 5 is the strongest.**

5) Place treatment window on treated area flatly, it is ready to press lighting button when light emitting indicators on. During operation, make sure no shots overlapping or missing on treated area, repeat the whole treatment for 2-3 times.

6) Press turn on/off button to turn off device, then cut off the power.

7) Use cotton swab or cleaning cloth to clean and disinfect light exit window.

## PREPARING for USE

1) Before each treatment, clean treatment area and shave areas you intend to treat as long as hairs continue to grow back.

2) Find your skin type with Skin Tone Test System and see suitable energy level on Treatment Level and Times Comparison Table.

3) Wear glasses.

4) Turn on the device and begin to shot to treatment area.

5) If you need, you can apply cooling product.

6) Apply moisturizing or sun block products after treatment.

## SUGGESTIONS for POWER LEVEL SELECTING AND OPERATION TIMES

- Two successive shots on underarms area or people with light skin tone,
- When the power level is higher, the operation times should be cut

down properly.

- Slightly redness or evenly warm feeling may be felt on treated area which is a kind of normal phenomenon.

## TREATMENT PROTOCOL

- On average, 6-10 treatments,
- For the first 3 to 5 sessions, spaced every 2 to 3 weeks,
- For the 6-10 sessions, spaced every 4 to 6 weeks,
- After that treatment should be done if hairs have grown back, until the desired results are achieved,
- The treatment interval is different from each other, the next treatment can be taken after 30% hairs grow back.

### Treatment Level and Times Comparison Table

| SKIN TYPE                    | WHITE                         | LIGHT BEIGE                   | BEIGE                         | LIGHT BROWN                   | MODERATE BROWN                | DARK BROWN & BLACK |
|------------------------------|-------------------------------|-------------------------------|-------------------------------|-------------------------------|-------------------------------|--------------------|
| Skin Tone Test System Result | I                             | II                            | III                           | IV                            | V                             | VI                 |
| Under arms                   | Level 3-5<br>Repeat 3-4 times | Level 3-5<br>Repeat 3-4 times | Level 3-4<br>Repeat 2-3 times | Level 3-4<br>Repeat 2-3 times | Level 2-3<br>Repeat 2-3 times | Not Suitable       |
|                              | Level 3-5<br>Repeat 3-4 times | Level 3-5<br>Repeat 3-4 times | Level 2-4<br>Repeat 2-3 times | Level 2-4<br>Repeat 2-3 times | Level 1-3<br>Repeat 2-3 times |                    |
| Arms                         | Level 3-5<br>Repeat 3-4 times | Level 3-5<br>Repeat 3-4 times | Level 2-4<br>Repeat 2-3 times | Level 2-4<br>Repeat 2-3 times | Level 1-3<br>Repeat 2-3 times | Not Suitable       |
|                              | Level 3-5<br>Repeat 3-4 times | Level 3-5<br>Repeat 3-4 times | Level 2-4<br>Repeat 2-3 times | Level 2-4<br>Repeat 2-3 times | Level 1-3<br>Repeat 2-3 times |                    |
| Legs                         | Level 3-5<br>Repeat 3-4 times | Level 3-5<br>Repeat 3-4 times | Level 2-4<br>Repeat 2-3 times | Level 2-4<br>Repeat 2-3 times | Level 1-3<br>Repeat 2-3 times | Not Suitable       |
|                              | Level 3-4<br>Repeat 3-4 times | Level 3-4<br>Repeat 3-4 times | Level 2-3<br>Repeat 2-3 times | Level 2-3<br>Repeat 2-3 times | Level 1-2<br>Repeat 2-3 times |                    |
| Lip                          | Level 3-4<br>Repeat 3-4 times | Level 3-4<br>Repeat 3-4 times | Level 2-3<br>Repeat 2-3 times | Level 2-3<br>Repeat 2-3 times | Level 1-2<br>Repeat 2-3 times | Not Suitable       |
|                              | Level 3-4<br>Repeat 3-4 times | Level 3-4<br>Repeat 3-4 times | Level 2-3<br>Repeat 2-3 times | Level 2-3<br>Repeat 2-3 times | Level 1-2<br>Repeat 2-3 times |                    |
| Facial                       | Level 3-4<br>Repeat 3-4 times | Level 3-4<br>Repeat 3-4 times | Level 2-3<br>Repeat 2-3 times | Level 2-3<br>Repeat 2-3 times | Level 1-2<br>Repeat 2-3 times | Not Suitable       |
|                              | Level 3-4<br>Repeat 3-4 times | Level 3-4<br>Repeat 3-4 times | Level 2-3<br>Repeat 2-3 times | Level 2-3<br>Repeat 2-3 times | Level 1-2<br>Repeat 2-3 times |                    |

\*Energy levels and repeat times are suggested for each treatment.

NOTE: First treat all treatment area without overlapping. Then, repeat

treatment as suggested on table above.

## CLEANING and MAINTANENCE

Clean device after each treatment for a long-lasting high performance.

Never wash device and any part of it.

Never use scouring pads, abrasive cleaning agents or aggressive liquids such as petrol or acetone to clean skin contact surface.

Switch off appliance and unplug it before cleaning.


If necessary, clean outside of the appliance with a soft cleaning cloth.

Moisten a soft lint-free cloth with a few drops of water and clean the skin contact surface and the light exit window.

## SIMPLE TROUBLE-SHOOTING

When fault occurs during the usage of the instrument, users can follow the following guidance to eliminate some simple faults. If the problem can not be solved after all steps are finished, please contact the Customer Service Department so as to solve faults as soon as possible.

- Failure of start-up: Please check the adapter whether connected well.
- Failure of lighting:
  - Please confirm whether the temperature of device is too high, if temperature is too high, the over-protection will start.
  - Please make sure the treatment window touches skin properly and whether the emitting light indicator turns light.
- Failure of energy level adjustment: Please check whether adjust energy levels more repeated more repeatedly, if so, please operate device after three minutes.
- Light emitting indicator on the two sides is always on: Please check whether the induction part of the head of the instrument is placed in the area in where the metal parts or induction capacitance effects exist.
- In case of abnormal situation or failure of fault elimination: Please contact our company or our local office and inform the purchase date. We will give you a reply as soon as possible.
  - Device doesnt pulse after a while: Over temperature protection systems automatically becomes activated in case of overheating and device stops pulsing for a while. Please wait until device becomes cool down.

Product Life 7 years

### Technical Specifications

12V DC 3A

Adaptor Input: 100-240V AC 50/60Hz 2A Max.

Adaptor Output: 12V DC 3A

**The product is in compliance with the EEE Directives.**

## WEEE

This product contains recyclable materials that compliance with the WEEE directives. Do not dispose this product as unsorted municipal waste. Please contact your local municipality for the nearest collection point.

**Manufacturer:**

**SHENZHEN GSD Tech CO., LTD.**

Building A, JUNSD Hi-Tech Park, West of Bao'an Rd. Watch & Clock

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